



Sartopore® 2 0.45 µm

Bioburden & Particle Reductive Filter Cartridges



Description

Sartopore® 2 0.45 µm rated filter cartridges are designed for bioburden reduction and particle removal from a broad range of pharmaceutical products. They offer extremely high flow rates and total throughputs and are therefore ideally suited for membrane prefiltration of aqueous solutions and highly viscous, difficult to filter pharmaceutical products.

Applications

Typical applications include bioburden reduction and particle removal from:

- Buffers
- Biological Fluids
- Ophthalmics
- LVP
- Antibiotics
- Bulk pharmaceutical products

Compatibility

Featuring a unique hydrophilic polyethersulfone membrane, Sartopore® 2 0.45 µm cartridges are compatible with solutions from pH 1 to pH 14 and are unaffected by numerous steam sterilization cycles. They are ideally suited for filtration of solutions with high | low pH and for multiple SIP | CIP cycles.

Performance

Sartopore® 2 0.45 µm cartridges provide exceptional high flow rates, resulting in economical sizing of filtration systems. Due to the "built-in prefiltration" by a 0.8 µm membrane, Sartopore® 2 0.45 µm rated cartridges offer outstanding total throughputs.

Wettability

Sartopore® 2 cartridges can be easily wetted out for integrity testing even after drying at 80 °C for 12 hours.

Microbiological Retention

Sartopore® 2 0.45 µm rated filter cartridges are validated for removal of *Serratia marcescens* with a Log Reduction Value (LRV) of 7 according to HIMA and ASTM F-838-05 guidelines.

Quality Control

Each individual element is integrity tested by Diffusion and Bubble Point test prior to release, assuring absolute reliability.

Documentation

Sartopore® 2 cartridges are designed, developed and manufactured in accordance with an ISO 9001 certified Quality Management System. A Validation Guide and Extractables Guide are available for compliance with regulatory requirements.

Specifications

Materials

Prefilter Membrane	Polyethersulfone, asymmetric
Endfilter Membrane	Polyethersulfone, asymmetric
Support Fleece	Polypropylene
Core	Polypropylene
End Caps	Polypropylene
O-Rings	Silicone (optional EPDM or Viton)

Pore Size Combination

0.8 µm + 0.45 µm

Available Sizes | Filtration Area

Cartridges

Size 0	5"	0.3 m ² 3.2 ft ²
Size 1	10"	0.6 m ² 6.5 ft ²
Size 2	20"	1.2 m ² 12.9 ft ²
Size 3	30"	1.8 m ² 19.4 ft ²

Mini Cartridges

Size 7	0.05 m ² 0.5 ft ²
Size 8	0.1 m ² 1.1 ft ²
Size 9	0.2 m ² 2.2 ft ²

Available Adapters Cartridges

- 21: double open end with flat gaskets
- 25: 2 Flange Bayonet adapter with 226 double o-ring
- 27: Bayonet adapter with 222 double o-ring
- 28: 3 Flange Bayonet adapter with 222 double o-ring

Available Adapters Mini Cartridges

- 15: Bayonet adapter with o-ring
- 18: Plug adapter with double o-ring

Operating Parameters

Max. Allowable Differential Pressure	5 bar 72.5 psi at 20°C
	2 bar 29 psi at 80°C
Max. Allowable Back Pressure	2 bar 29 psi at 20°C

Specifications

Extractables

Sartopore® 2 0.45 µm rated filter cartridges meet, or exceed the requirements for WFI quality standards set by the current USP.

Regulatory Compliance

Individually integrity tested

Integrity test correlated to HIMA/ASTM F 838-05 Bacteria Challenge Test using *Serratia marcescens*

Non-pyrogenic according to USP Bacterial Endotoxins

Meets USP Plastics Class VI biological reactivity test, in vivo

Non-fiber releasing according to 21 CFR

Sterilization

In-Line Steam Sterilization:

134 °C, 20 min. at max differential pressure of 0.5 bar | 7.25 psi

Autoclaving:

134 °C, 2 bar | 29 psi, 30 min

Sterilization Cycles

In-Line Sterilization: Min. 25
Autoclaving: Min. 25

Technical References

Validation Guide:

SPK 5732-e

Extractables Guide:


SPK5731-e

Ordering Information

Order Code	Pore size [µm]	Pack size [Pieces]	Test Pressure [bar psi]	Max. Diffusion [ml/min]	Min. Bubble Point [bar psi]
Standard Cartridges					
544**06G0	0.45	1	1.7 25	7	2.2 32
544**06G1	0.45	1	1.7 25	12	2.2 32
544**06G2	0.45	1	1.7 25	24	2.2 32
544**06G3	0.45	1	1.7 25	36	2.2 32
Mini Cartridges					
544**06G7-----B	0.45	5	1.7 25	3	2.2 32
544**06G8-----B	0.45	5	1.7 25	4	2.2 32
544**06G9-----B	0.45	5	1.7 25	6	2.2 32


** : Adapter Type

Water Flow Rates for Standard Cartridges

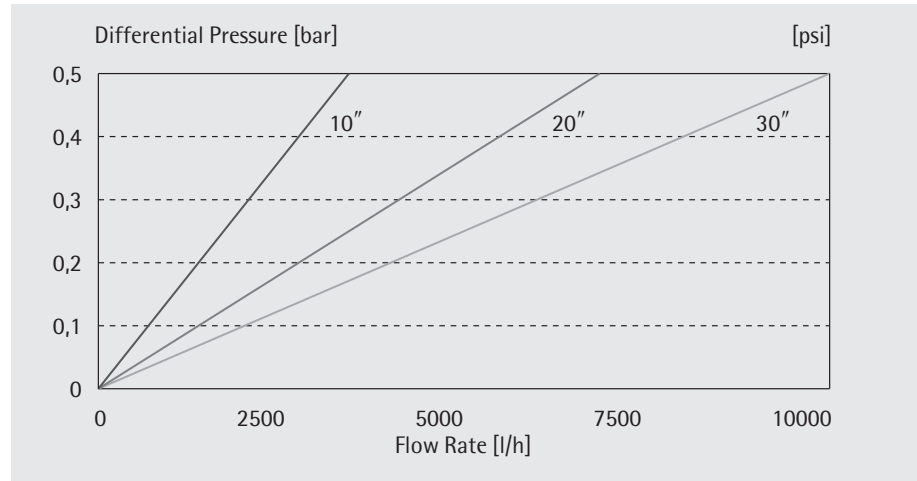

Standardized at 20 °C

Total Throughput Comparison

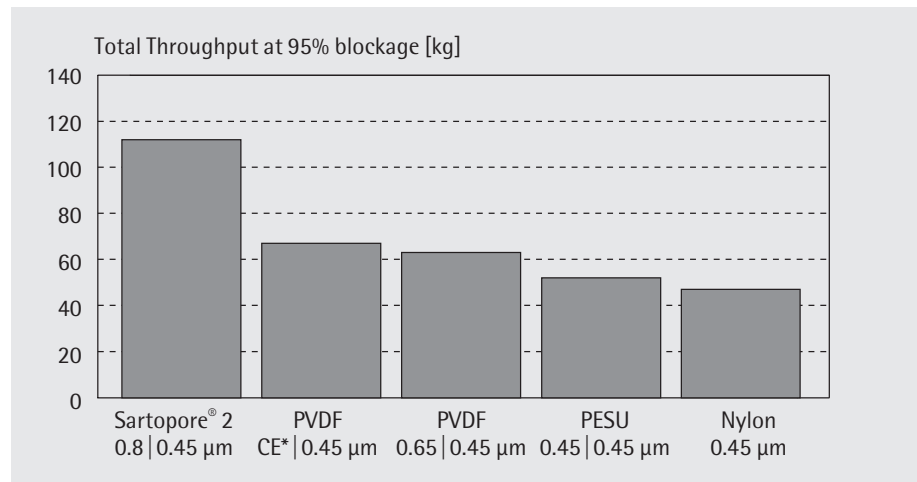

10" Cartridges

* Cellulose Ester prefilter

Sartorius Stedim Biotech GmbH
August-Spindler-Strasse 11
37079 Goettingen, Germany

Phone +49.551.308.0
Fax +49.551.308.3289
www.sartorius-stedim.com

USA Toll-Free +1.800.368.7178
UK +44.1372.737159
France +33.442.845600
Italy +39.055.63.40.41
Spain +34.90.2110935
Russian Federation +7.812.327.5.327
Japan +81.3.5639.9981

Specifications subject to change without notice. Printed and copyrighted by Sartorius Stedim Biotech GmbH. | W
Publication No.: SPK2039-e13057
Order No.: 85030-513-44
Ver. 05 | 2013